



**IBA Lifesciences GmbH**  
Rudolf-Wissell-Str. 28  
37079 Goettingen  
Germany  
Tel.: +49 (0) 551-5 06 72-0  
E-mail: [info@iba-lifesciences.com](mailto:info@iba-lifesciences.com)  
[www.iba-lifesciences.com](http://www.iba-lifesciences.com)

## Data Sheet

### pESG-IBA143

Cat. No.: 5-4543-001

Version: 3.0  
Revision Date: 27.07.2021

<b>Description</b>	<b>StarGate Acceptor Vector for stable and non-replicative transient expression in mammalian cells. It carries the human cytomegalovirus (CMV) immediate-early promoter for high-level expression and the SV40 ori for episomal replication in cell lines that are latently infected with SV40 or that express the SV40 large T antigen (e.g., HEK293T, COS-1, COS-7). The expressed recombinant protein will be localized in the cytoplasm.</b>
<b>Affinity tag</b>	Twin-Strep-tag <sup>®</sup> is fused to the C-terminus and 6xHistidine-tag is fused to the N-terminus of the recombinant protein.
<b>Cloning Strategy</b>	Cloning into StarGate Acceptor Vectors has to be done with the restriction enzyme Esp3I. There is no Multiple Cloning Site (MCS) available that can be used for the integration of the gene of interest instead (see manual).
<b>Resistance</b>	Ampicillin: for selection of transformed <i>E. coli</i> cells Neomycin: for selection of stable cell lines
<b>Form</b>	5 µg, dissolved in 20 µl TE buffer, pH 8.0: 10 mM Tris/HCl, 1 mM EDTA
<b>Concentration</b>	250 ng/µl
<b>Stability</b>	12 months after shipping
<b>Storage</b>	recommended: 2-8 °C for frequent usage, -20 °C for long-term storage
<b>Shipping</b>	room temperature
<b>Hazards</b>	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.

#### For research use only


##### Trademark information

The owners of trademarks marked by “®” or “TM” are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

##### Important licensing information


This product is covered by intellectual property (IP) rights and on completion of the sale IBA Lifesciences grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at [info@iba-lifesciences.com](mailto:info@iba-lifesciences.com) or at IBA Lifesciences GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.

## Expression cassette of pESG-IBA143

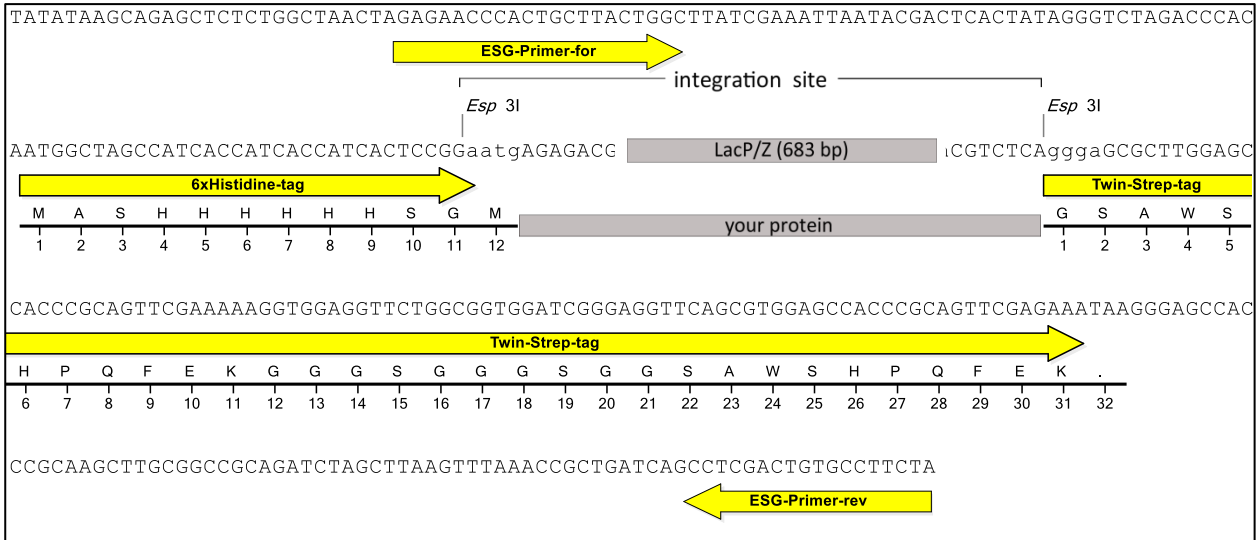

LacP/Z cassette = contains LacZ alpha fragment under control of a separate promoter, which allows alpha complementation of *LacZ* mutations such as *LacZ* $\Delta$ M15 as in *E. coli* DH5 $\alpha$  or TOP10.

your protein = after StarGate cloning using *Esp31* your gene of interest will be located here

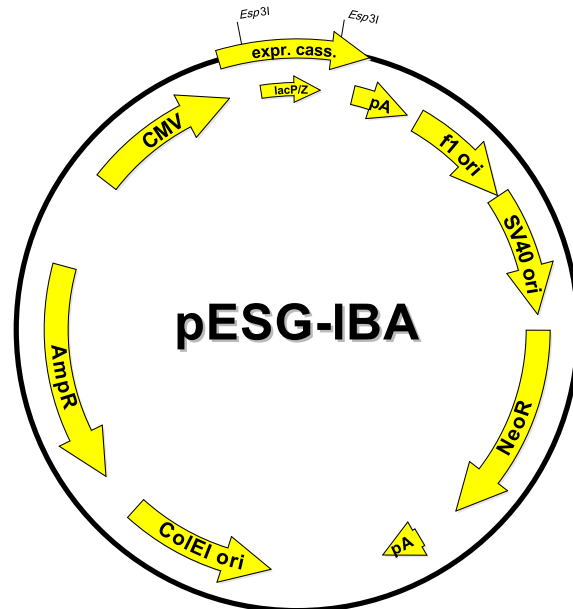

Features	from bp	to bp	Sequencing primer
f1 origin	259	687	ESG-Primer-for
SV40 ori	692	1035	
Neomycin resistance gene	1097	1891	5' - GAGAACCCACTGCTTACTGGC -3'
ColEI ori	2637	3222	ESG-Primer-rev
Ampicillin resistance gene	3393	4253	
CMV promoter	4621	5208	5' - TAGAAGGCACAGTCGAGG -3'
forward primer binding site	5221	5241	
6xHistidine-tag	5284	5316	
LacZ alpha fragment	5545	5946	
Twin-Strep-tag®	6010	6102	
reverse primer binding site	6164	6181	
polyA signal sequence	1	213	
total vector length		6181	