

Data Sheet

pESG-IBA35

Cat. No.: 5-4435-001

Version: 2.2

Lot No.: 4435-

Revision Date: 04.03.2020

Description	StarGate® Acceptor Vector designed for high-level stable and non-replicative transient expression in most mammalian hosts containing the following elements: <ul style="list-style-type: none"> • Human cytomegalovirus (CMV) immediate-early promoter for high-level expression in a wide range of mammalian cells • Neomycin resistance gene for selection of stable cell lines • Episomal replication in cell lines that are latently infected with SV40 or that express the SV40 large T antigen (e.g. COS-1, COS-7) • The expressed recombinant protein will be localized in the cytoplasm.
Affinity tag	6xHistidine-tag for the purification of recombinant protein via Ni-NTA resins. The 6xHistidine-tag is fused to the N-terminus of the recombinant protein.
Resistance	Ampicillin
Form	5 µg, dissolved in 20 µl TE buffer, pH 8,0: 10 mM Tris-HCl, 1 mM EDTA
Concentration	250 ng/µl
Stability	12 months after shipping
Storage	recommended: 2-8 °C for frequent usage, -20 °C for long-term storage
Shipping	room temperature
Hazards	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.

Note: The sequences have been compiled from information in the sequence database, published literature, and other sources, together with partial sequences obtained by IBA, however, the vectors have not been completely sequenced.



Go digital and help the environment. Please download all up-to-date manuals, protocols and other material from <http://www.iba-lifesciences.com>.

For research use only


Important licensing information

This product is based on StarGate, 6xHistidine and CMV promoter technologies covered by intellectual property (IP) rights and on completion of the sale IBA grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at info@iba-lifesciences.com or at IBA GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.

Trademark information


The owners of trademarks marked by “®” or “TM” are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

Expression cassette of pESG-IBA35

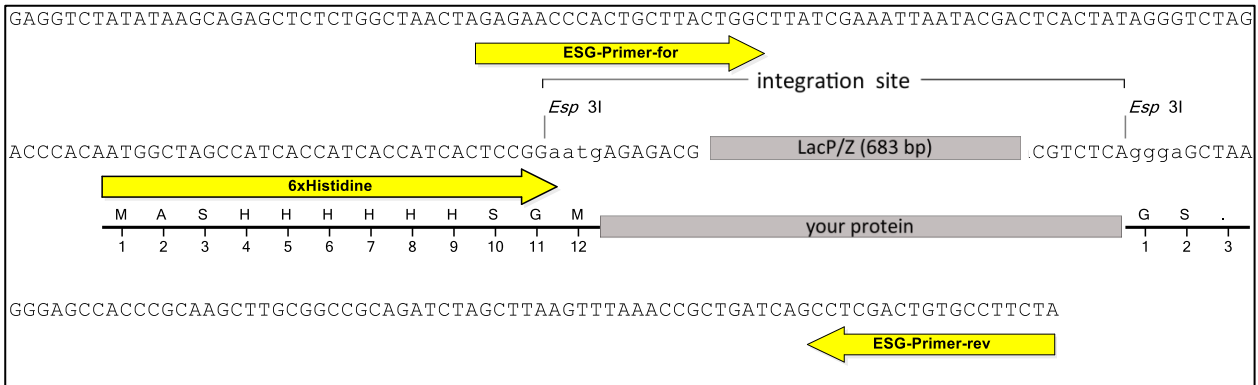

LacP/Z cassette = contains LacZ alpha fragment under control of a separate promoter, which allows alpha complementation of *LacZ* mutations such as *LacZΔM15* as in *E. coli* DH5α or TOP10.

your protein = after StarGate cloning using *Esp3I* your gene of interest will be located here

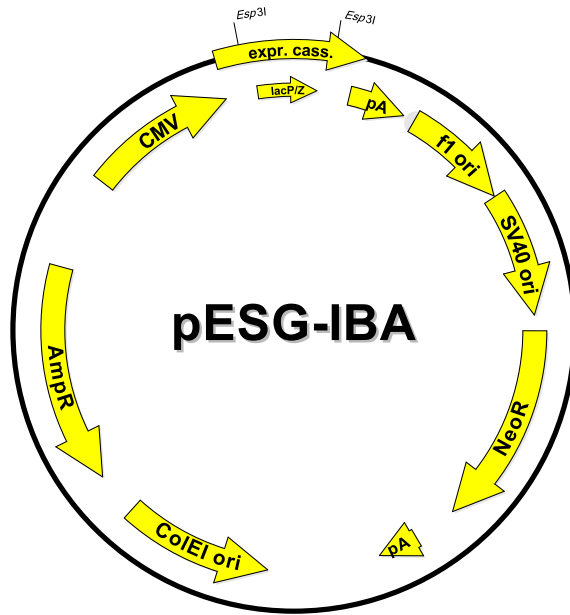

Features	from bp	to bp	Sequencing primer
f1 origin	259	687	ESG-Primer-for (Cat. No. 5-0000-121)
SV40 ori	692	1035	
Neomycin resistance gene	1097	1891	5' - GAGAACCCACTGCTTACTGGC -3'
ColEI ori	2637	3222	ESG-Primer-rev (Cat. No. 5-0000-122)
Ampicillin resistance gene	3393	4253	
CMV promoter	4621	5208	5' - TAGAAGGCACAGTCGAGG -3'
forward primer binding site	5221	5241	
6xHistidine-tag	5284	5316	
LacZ alpha fragment	5545	5946	
reverse primer binding site	6077	6094	
polyA signal sequence	1	213	
total vector length		6094	