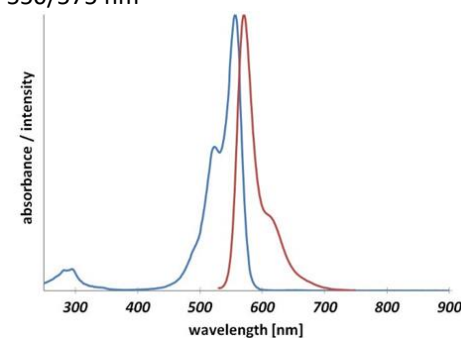


## Data Sheet

# StrepMAB-Classic DY-549

Cat. No.: 2-1566-050

Version: 2.3  
 Revision Date: 02.06.2022

<b>Description</b>	<b>StrepMAB-Classic DY-549 is a highly selective, murine monoclonal antibody for detection of Strep-tag®II and Twin-Strep-tag® fusion proteins in applications like western blots, immunohistochemistry, and immunocytochemistry. It is conjugated with the fluorescent dye DY-549.</b>
<b>Specificity</b>	Strep-tag®II, Twin-Strep-tag®
<b>Product Type</b>	Monoclonal antibody (mouse IgG1)
<b>Form</b>	Suspension in PBS; contains 0.01% NaN <sub>3</sub> , 1 mg/ml BSA
<b>Concentration</b>	0.5 mg/ml
<b>Amount</b>	100 µl
<b>Stability</b>	6 months after shipping (If the products contain a precipitate, we recommend centrifugation before use.)
<b>Storage</b>	2-8 °C, protect from light
<b>Shipping</b>	Room temperature
<b>Hazards</b>	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.
<b>Absorption/ Emission</b>	550/575 nm  Data kindly provided by Dyomics GmbH
<b>Application</b>	If StrepMAB-Immo DY-549 conjugate does not work, e.g., due to a missing SA-linker, use StrepMAB-Classic DY-549 conjugate in ELISA with a dilution between 1:100 to 1:1000.

### For research use only

#### Trademark information

The owners of trademarks marked by “®” or “TM” are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

#### Important licensing information

This product is covered by intellectual property (IP) rights and on completion of the sale IBA Lifesciences grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at [info@iba-lifesciences.com](mailto:info@iba-lifesciences.com) or at IBA Lifesciences GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.