



**BA Lifesciences GmbH**  
Rudolf-Wissell-Str. 28  
37079 Goettingen  
Germany  
Tel.: +49 (0) 551-5 06 72-0  
E-mail: [info@iba-lifesciences.com](mailto:info@iba-lifesciences.com)  
[www.iba-lifesciences.com](http://www.iba-lifesciences.com)

## Data Sheet

# pASG-IBA64

Cat. No.: 5-4064-001

Version: 3.0  
Revision Date: 26.07.2021

<b>Description</b>	<b>StarGate Acceptor Vector for bacterial expression. The expression cassette is under transcriptional control of the tetracycline promoter/operator. The expressed recombinant protein will be localized in the cytoplasm.</b>
<b>Affinity tag</b>	Strep-tag® II is fused to the N-terminus and FLAG-tag is fused to the C-terminus of the recombinant protein.
<b>Cloning Strategy</b>	Cloning into StarGate Acceptor Vectors has to be done with the restriction enzyme Esp3I. There is no Multiple Cloning Site (MCS) available that can be used for the integration of the gene of interest instead (see manual).
<b>Expression strain</b>	Any <i>E. coli</i> strain. The <i>tet</i> -promoter works independently from the genetic background of <i>E. coli</i> .
<b>Bacterial Expression</b>	Expression is induced upon addition of 200 µg anhydrotetracycline per 1 liter <i>E. coli</i> shaking culture ( $A_{550} = 0.5$ ).
<b>Resistance</b>	Ampicillin
<b>Form</b>	5 µg, dissolved in 20 µl TE buffer, pH 8.0: 10 mM Tris/HCl, 1 mM EDTA
<b>Concentration</b>	250 ng/µl
<b>Stability</b>	12 months after shipping
<b>Storage</b>	recommended: 2-8 °C for frequent usage, -20 °C for long-term storage
<b>Shipping</b>	room temperature
<b>Hazards</b>	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.

### For research use only


#### Trademark information

The owners of trademarks marked by “®” or “TM” are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

#### Important licensing information


This product is covered by intellectual property (IP) rights and on completion of the sale IBA Lifesciences grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at [info@iba-lifesciences.com](mailto:info@iba-lifesciences.com) or at IBA Lifesciences GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.

## Expression cassette of pASG-IBA62

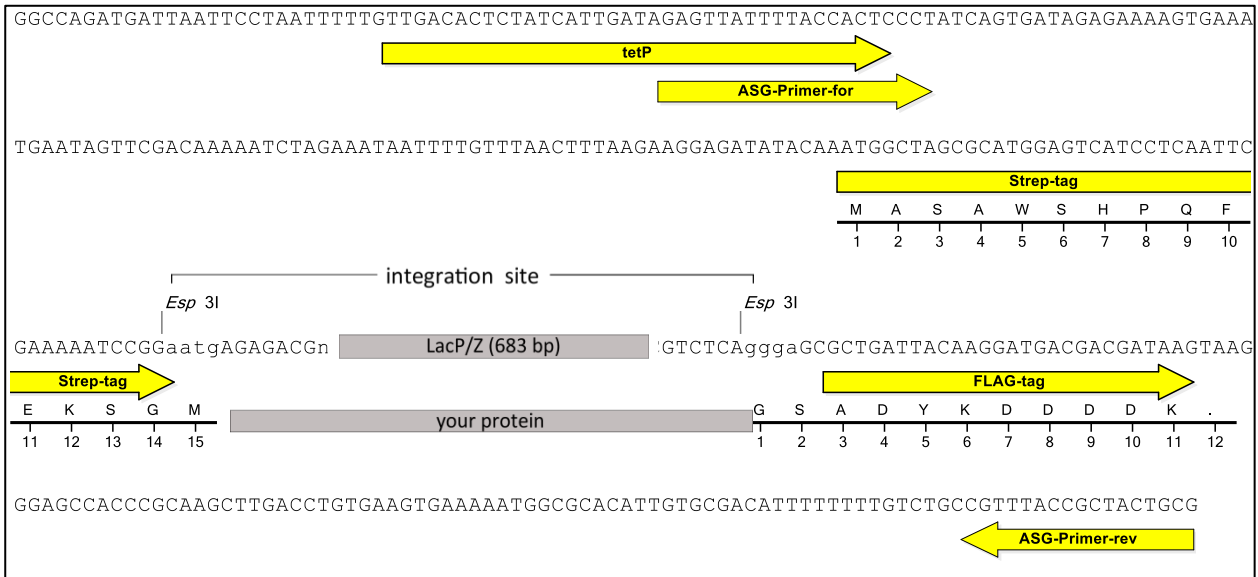

LacP/Z cassette = contains LacZ alpha fragment under control of a separate promoter, which allows alpha complementation of *LacZ* mutations such as *LacZ*ΔM15 as in *E. coli* DH5α or TOP10.

your protein = after StarGate cloning using *Esp3I* your gene of interest will be located here

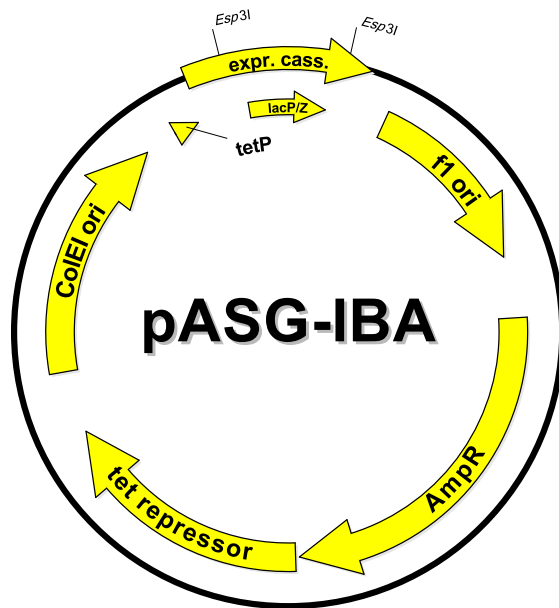

Features	from bp	to bp	Sequencing primer
f1 origin	13	451	<b>ASG-Primer-for</b>
AmpR resistance gene	600	1460	
Tet-repressor	1470	2093	5' - GAGTTATTTTACCACTCCCT - 3'
ColEI ori	2246	2834	<b>ASG-Primer-rev</b>
Tet promoter	2939	2975	
forward primer binding site	2959	2978	5' - CGCAGTAGCGGTAAACG - 3'
Strep-tag®II	3062	3103	
LacZ alpha fragment	3332	3733	
FLAG-tag	3803	3829	
reverse primer binding site	3903	3919	
total vector length		3919	