



IBA Lifesciences GmbH
Rudolf-Wissell-Str. 28
37079 Goettingen
Germany
Tel.: +49 (0) 551-5 06 72-0
E-mail: info@iba-lifesciences.com
www.iba-lifesciences.com

Data Sheet

pLSG-IBA164

Cat. No.: 5-4964-001

Version: 3.0
Revision Date: 27.07.2021

Description	StarGate Acceptor Vector for high-level expression of target proteins in insect cells. The gene transfer into the polyhedrin gene locus of AcMNPV DNA is achieved by homologous recombination and the vector carries a polyhedrin promoter. Co-transfection with BacPAK6 linearized AcMNPV DNA (Clontech) or with circular flashBAC modified AcMNPV DNA (Oxford Expression Technologies) allows the generation of recombinant baculovirus at very high efficiency through reconstitution of an essential gene (ORF 1629) and elimination of wild-type virus to great extent. The expressed recombinant protein will be localized in the cytoplasm.
Affinity tag	Twin-Strep-tag® is fused to the N-terminus and FLAG-tag is fused to the C-terminus of the recombinant protein.
Cloning Strategy	Cloning into StarGate Acceptor Vectors has to be done with the restriction enzyme Esp3I. There is no Multiple Cloning Site (MCS) available that can be used for the integration of the gene of interest instead (see manual).
Resistance	Ampicillin: for selection of transformed <i>E. coli</i> cells
Form	5 µg, dissolved in 20 µl TE buffer, pH 8.0: 10 mM Tris/HCl, 1 mM EDTA
Concentration	250 ng/µl
Stability	12 months after shipping
Storage	recommended: 2-8 °C for frequent usage, -20 °C for long-term storage
Shipping	room temperature
Hazards	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.

For research use only


Trademark information

The owners of trademarks marked by “®” or “TM” are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

Important licensing information


This product is covered by intellectual property (IP) rights and on completion of the sale IBA Lifesciences grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at info@iba-lifesciences.com or at IBA Lifesciences GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.

Expression cassette of pLSG-IBA164

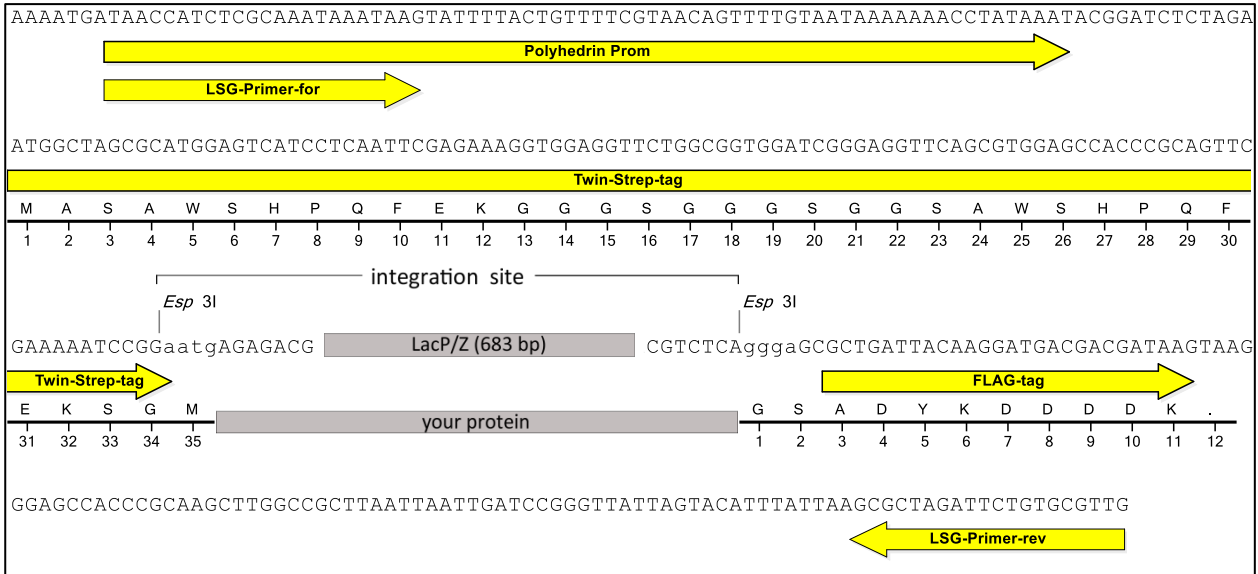

LacP/Z cassette = contains LacZ alpha fragment under control of a separate promoter, which allows alpha complementation of *LacZ* mutations such as *LacZΔM15* as in *E. coli* DH5α or TOP10.

your protein = after StarGate cloning using *Esp3I* your gene of interest will be located here

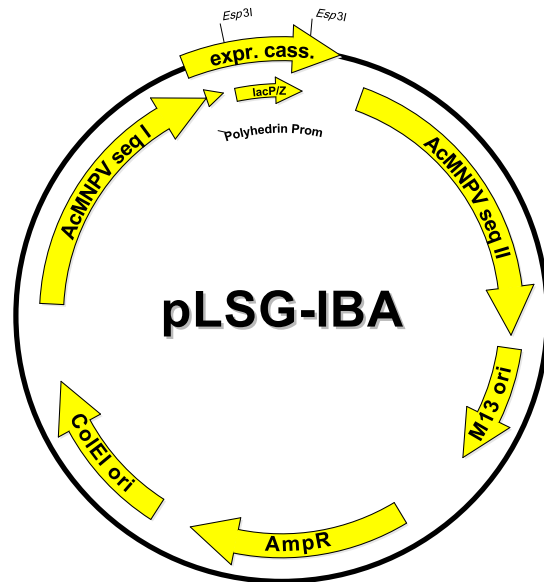

Features	from bp	to bp	Sequencing primer
AcMNPVseq II	1	1395	LSG-Primer-for
M13 ori	1447	1920	
Ampicillin resistance gene	2251	3111	5' - TAACCATCTCGCAAATAAATAAG -3' LSG-Primer-rev 5' - CAACGCACAGAATCTAGCGC -3'
ColEI ori	3259	3902	
AcMNPVseq I	4211	5357	
Polyhedrin promoter	5286	5355	
forward primer binding site	5286	5308	
Twin-Strep-tag®	5369	5470	
LacZ alpha fragment	5699	6100	
FLAG-tag	6170	6196	
reverse primer binding site	6262	6281	
total vector length		6281	